


HybridKnife®



High pressure tissue-selective hydrodissection

The **01 HybridKnife®** combines the **02 VIO®** electrosurgical system and the **03 ERBEJET® 2** hydrosurgery system in a single hybrid instrument for endoscopic procedures.

- Needle-free injection on demand
- Cutting and coagulation with proprietary VIO® technology

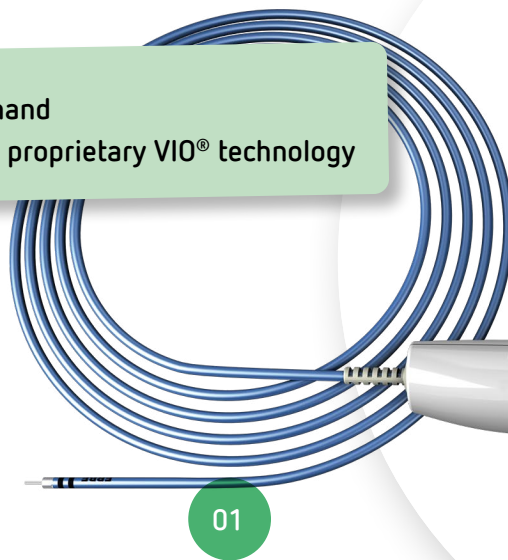

INJECT - CUT - COAGULATE: NO NEED FOR INSTRUMENT EXCHANGE

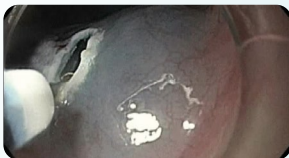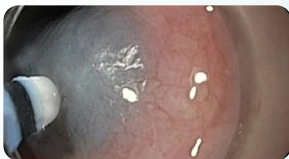


- **High pressure tissue-selective hydrodissection** defines tissue planes and creates a protective saline cushion to further dissipate energy and mitigate unintended thermal damage
- **Consistent and reliable protective cushion**, even in challenging anatomy, e.g. submucosal fibrosis, tortuous folds

- **Injection on demand** can be utilized to isolate and expose vasculature
- **Two distal tip designs** with a common 0.7 mm shaft geometry facilitates enhanced coagulation

HybridKnife®

Product data

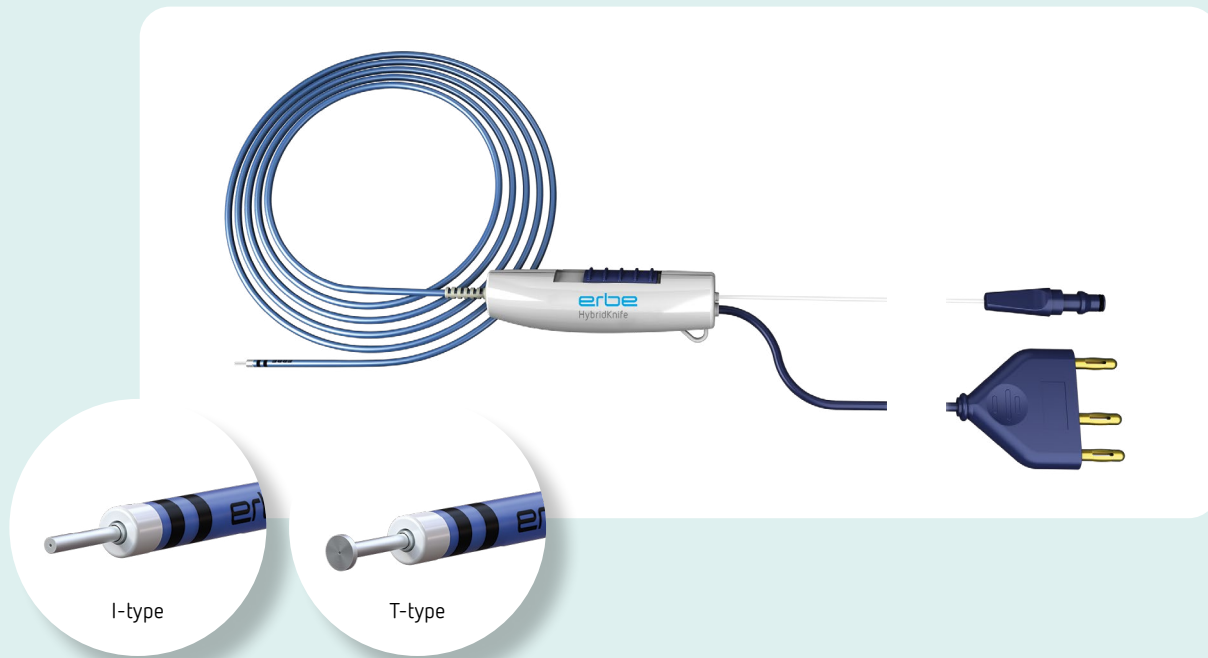

Single-use instruments and accessories

REF No.	
20150-260	HybridKnife®, T-type, I-jet O.D. 2.3 mm, length 1.9 m
20150-261	HybridKnife®, I-type, I-jet O.D. 2.3 mm, length 1.9 m
20150-220	ERBEJET® 2 probe, flexible, OD 1.3 mm, length 2.2 m (box of 5)

Units, modules and accessories

REF No.	
10150-000	ERBEJET® 2 hydrosurgical unit
10160-000	VIO® 3 electrosurgical unit
20150-300	ERBEJET® 2 pump cartridge
20150-102	ERBEJET® 2 one-pedal footswitch w/ ReMode®, AP & IP X8 equipment
20180-000	VIO® cart system carrier for VIO®
20180-010	VIO® cart wire basket
20180-040	Bracket/side rail, length 260 mm
20180-042	Bracket/side rail, length 390 mm
20180-141	Fastening set ERBEJET® 2/APC® 3 to VIO® cart
20180-144	Fastening set VIO® 3 and APC® 3 to ERBEJET® 2

Erbe USA Incorporated
2225 Northwest Parkway
Marietta, GA 30067
USA

Phone 800 778 3723
Fax 770 955 2577
sales@erbe-usa.com
erbe-usa.com